

# Product information



**Win 102 BTE**

**Win 112 BTE**


**Win 105 BTE DM**



Win hearing instruments have the ideal combination of affordable digital technology, superior product quality, and easy to use fitting software.

### Features to provide client satisfaction

- 3 channels
- Adaptive Noise Reduction
- Directional microphone system with active microphone matching (on 105 BTE DM)
- Soft Noise Management™
- Feedback Management
- Programmable beep indicators
- Telecoil boost
- Push button control of hearing instrument modes
- Rotary volume control
- FM compatibility
- ON/OFF function on battery drawer

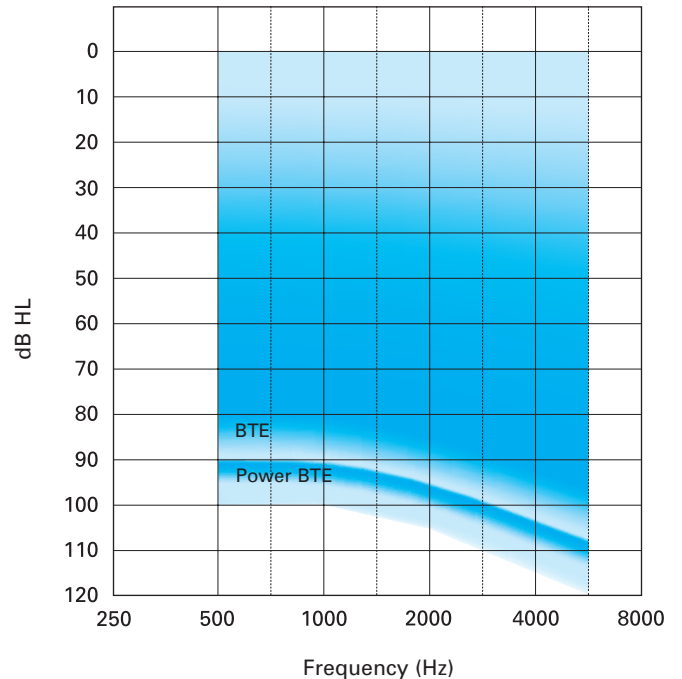

### **Win 102 BTE**

A compact BTE with a size 13 battery, for mild to severe hearing losses.



### **Win 105 BTE DM**

A compact BTE with a size 13 battery and directional microphone technology, for mild to severe hearing losses.



### **Win 112 BTE**

A compact power BTE with a size 13 battery, for mild to severe/profound hearing losses.

# Win information

Win BTE								
			102 BTE		112 BTE		105 BTE DM	
			IEC 60118-7 (2cc)	IEC 60118-0 (Ears.)	IEC 60118-7 (2cc)	IEC 60118-0 (Ears.)	IEC 60118-7 (2cc)	IEC 60118-0 (Ears.)
Output	OSPL 90, Peak	dB SPL	123	128	132	136*	123	127
	OSPL 90, 1600 Hz,	dB SPL	115	123	122	130	113	121
	HFA OSPL 90, ANSI	dB SPL	118	—	126	—	113	—
Gain	Full-On Gain, Peak	dB	60	64	71	76	60	66
	Full-On Gain, 1600 Hz	dB	51	59	60	68	49	59
	HFA Full-On Gain, ANSI	dB	55	—	65	—	50	—
	Reference Test Gain, IEC	dB	39	48	47	55	38	46
	Reference Test Gain, ANSI	dB	41	—	49	—	36	—
Current	Quiescent Current	mA	1.0	1.0	1.1	1.1	1.0	1.0
	Operating Current, IEC	mA	1.0	1.0	1.1	1.1	1.0	1.0
	Operating Current, ANSI	mA	1.0	—	1.2	—	1.0	—
	Battery Type		13		13		13	
Distortion	500/800/1600 Hz typ., IEC	%	<1 / <1 / <1	<1 / <1 / <1	<1 / <1 / <1	<1 / <1 / <1	<1 / <1 / <1	<2 / <1 / <1
	500/800/1600 Hz typ., ANSI	%	<1 / <1 / <1	—	<1 / <1 / <1	—	<1 / <1 / <1	—
General Information	Frequency Range, ANSI	Hz	130–6250	—	140–6000	—	100–7000	—
	Equiv. Input Noise, IEC/ANSI <sup>1)</sup>	dB	19/14		20/15		19/19	
	Telecoil 1 mA/m 1600 Hz, IEC	dB SPL	85	92	88	96	82	91
	Telecoil HFA SPLITS, ANSI	dB SPL	103	—	111	—	95	—
Additional Information	Colors		beige		beige		beige	
	Earhook		8 dB (std.), 0 dB (opt)		8 dB (std.), 0 dB (opt)		8 dB (std.), 0 dB (opt)	
	Child's Earhook		8 dB (opt), 0 dB (opt)		8 dB (opt), 0 dB (opt)		8 dB (opt), 0 dB (opt)	
	Telecoil Boost		yes		yes		yes	
	Push Button Modes		M, T or MT / Mute		M, T or MT / Mute		M, DM, T or MT / Mute	
	FM Communication		available		available		available	
	DAI		yes		yes		yes	
	Volume Control		yes		yes		yes	
	AI-DI	0°/max. dB	—		—		2.9/3.1	

<sup>1)</sup> Technical data measured with expansion, corresponding to Soft Noise Management™ level 3.

All measurements are made according to IEC 60118 if not otherwise mentioned. ANSI refers to ANSI S3.22-1996.

The Full-On Gain setting can be programmed into the instrument from OASIS plus for verification purposes.

All measurements are based on earhook without filter.

Directionality measured on Kemar®. AI-DI relative to 0° and relative to maximal polar sensitivity, calculated for the horizontal plane.

\* Special care should be exercised in selecting and fitting a hearing aid whose maximum sound pressure level exceeds 132 dB, as measured with an IEC 60711: 1981 occluded ear simulator, because there may be a risk of impairing the remaining hearing of the hearing instrument user.

## Programming

Win is programmed with OASIS plus, Bernafon's outstanding fitting software.

OASIS plus is a NOAH compatible, MS-Windows based PC-Fitting software.

The connection is made through a HiPro or NOAHlink interface, using the appropriate programming cables and adapters.

System requirements:


NOAH 3 with HiPro or NOAHlink,

programming cables and programming adapters.



# Frequency responses


## IEC 60118-7 2cc Coupler (IEC 60126)

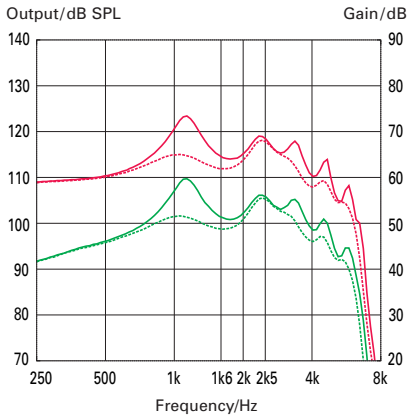

### Win 102 BTE



Output OSPL 90  
Full-On Gain

## IEC 60118-0 Earsimulator (IEC 60711)

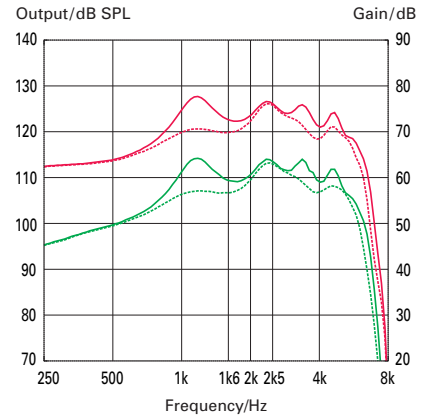



### Win 112 BTE

Output OSPL 90  
Full-On Gain

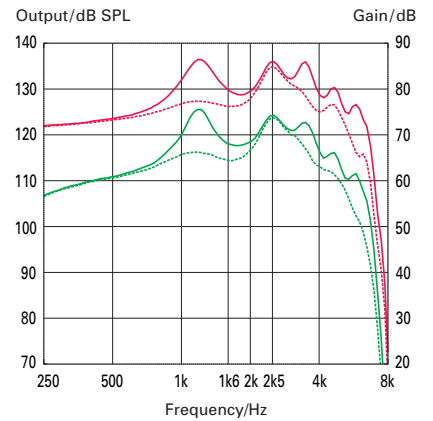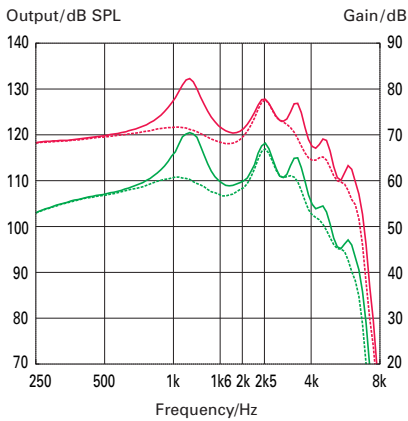

### Win 105 BTE DM

Output OSPL 90  
Full-On Gain

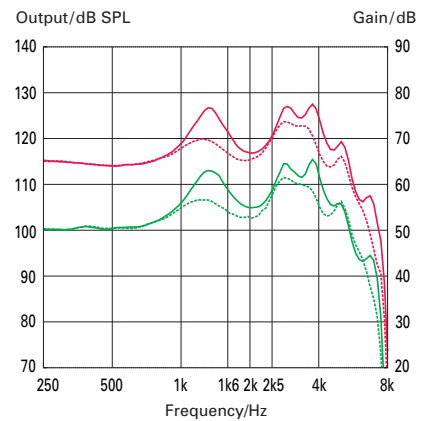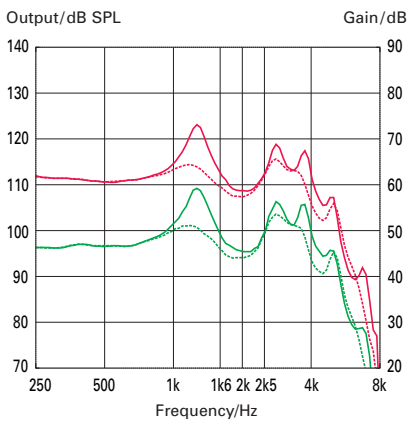

— Measurements based on earhook without filter  
- - - Measurements based on earhook with filter

# DAI/FM Communication



### FMA3 Adapter

Part No. 399-50-570-00

For LEXIS and other FM communication systems only



### DA13 Adapter

Part No. 399-50-580-00

For connection to Hi-Fi, TV set, PC, etc.


# Fitting in OASIS plus

*A high performance fitting software platform to ensure a perfect fit.*

All Win hearing instruments are fit with OASIS plus, the reliable hearing instrument fitting software from Bernafon.

Intelligent software makes your job easier.

- First fit in as few as 3 mouse clicks
- Fine control for each channel for greater adjustment precision
- 3 MPO controls
- Interactive fine-tuning to support the dialog with your client
- Efficient Feedback Management
- COSI for counseling and ensuring that your client's goals are met
- NOAHlink compatibility
- Standard cables and adapters
- Unified fitting flow with all Bernafon digital hearing instruments

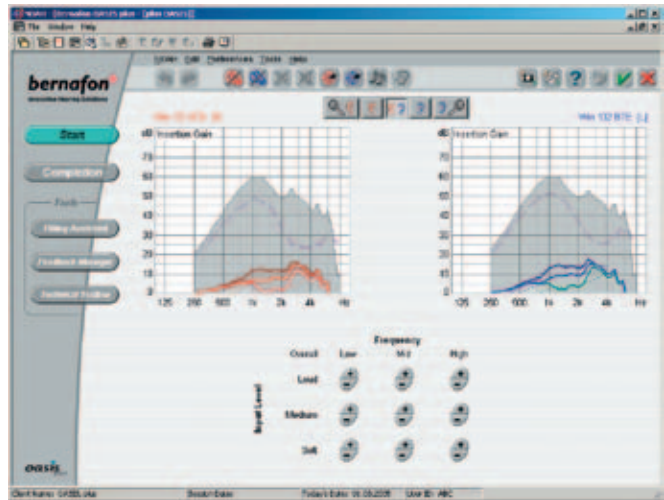

*Win and OASIS plus, the perfect combination of affordable hearing instruments and high performance fitting software.*

## Programming cable attachment



### **Programming cables:**

Nr. 2, NEW STANDARD (HiPro)  
Blue, left Ref. 384-20-033-00  
Red, right Ref. 384-20-032-00

### **Programming cables:**

Nr. 2, NEW STANDARD (NOAHlink)  
Blue, left Ref. 384-20-035-00  
Red, right Ref. 384-20-034-00



### **Programming adapter**

Programming adapter Ref. 390-01-151-00

### **World Headquarters**

#### **Switzerland**

Bernafon AG  
Morgenstrasse 131  
3018 Bern  
Phone +41 31 998 15 15  
Fax +41 31 998 15 90  
www.bernafon.com

#### **Australia**

Bernafon Australia Pty. Ltd.  
12/97 Castlemaine Street  
Milton QLD 4064  
Freecall 1800 809 111  
Phone +61 7 3250 0300  
Fax +61 7 3250 0372  
www.bernafon.com.au

#### **New Zealand**

Bernafon New Zealand Ltd.  
Level 1, Building F  
27-29 William Pickering Drive  
Albany, Auckland 0632  
Toll Free 0800 44 22 57  
Telephone +64 9 415 7917  
Fax +64 9 415 7916  
www.bernafon.co.nz

#### **United Kingdom**

Bernafon UK  
Cadzow Industrial Estate  
Off Low Waters Road  
Hamilton  
ML3 7QE Scotland  
Phone +44 1698 285 968  
Fax +44 1698 421 456  
www.bernafon.co.uk

**SWISS**   
Engineering

**bernafon**   
Your hearing • Our passion

# Product information



**Win 202 ITE**

**Win 322 ITC**

**Win 301 ITC/Win 302 ITC**

**Win 401 CIC/Win 411 MC**


**Win 315 ITC DM**



Win hearing instruments have the ideal combination of affordable digital technology, superior product quality, and easy to use fitting software.

## Features to provide client satisfaction

- 3 channels
- Adaptive Noise Reduction
- Directional microphone system with active microphone matching (on 315 ITC DM)
- Soft Noise Management™
- Feedback Management
- Programmable beep indicators
- Automatic telecoil (optional, except 401/411)
- Telecoil boost (optional, except 401/411)
- Push button control of hearing instrument modes

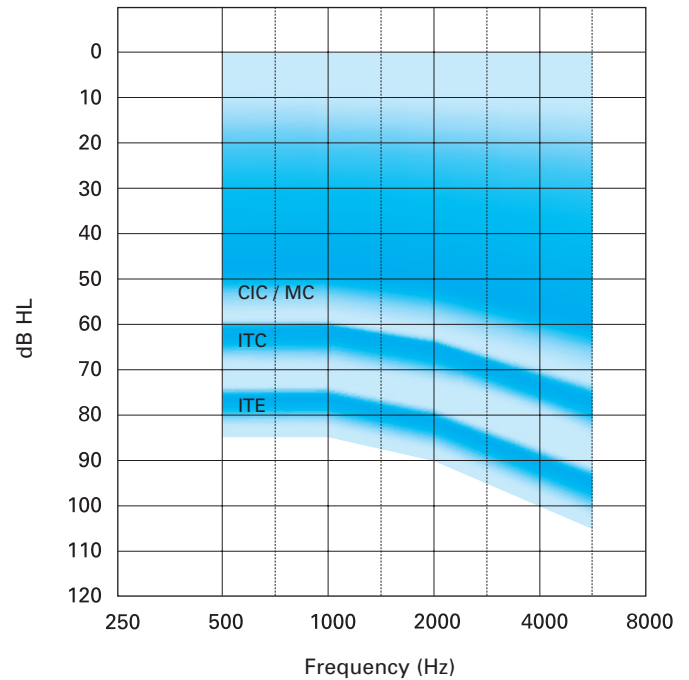

### Win 202 ITE

A power concha instrument with a size 13 battery and rotary volume control, for moderate to severe hearing losses.



### Win 322 ITC/315 ITC DM/302 ITC

A canal instrument with a size 312 battery and rotary volume control, for mild to moderately-severe hearing losses. Dual microphone technology on 315 ITC DM.



### Win 301 ITC

A small canal instrument with a size 312 battery, for mild to moderate hearing losses.



### Win 411 MC

Mini canal instrument with a size 10 battery, for mild to moderate hearing losses.



### Win 401 CIC

Completely in the canal instrument with a size 10 battery, for mild to moderate hearing losses.

# Win information

Win ITE								
		202 ITE		322 ITC		302 ITC		
		IEC 60118-7 (2cc)	IEC 60118-0 (Ears.)	IEC 60118-7 (2cc)	IEC 60118-0 (Ears.)	IEC 60118-7 (2cc)	IEC 60118-0 (Ears.)	
Output	OSPL 90, Peak	dB SPL	120	129	114	124	110	120
	OSPL 90, 1600 Hz	dB SPL	112	120	105	113	103	111
	HFA OSPL 90, ANSI	dB SPL	114	—	106	—	104	—
Gain	Full-On Gain, Peak	dB	53	62	44	56	40	50
	Full-On Gain, 1600 Hz	dB	43	52	33	41	31	39
	HFA Full-On Gain, ANSI	dB	45	—	34	—	32	—
	Reference Test Gain, IEC	dB	36	45	26	34	23	32
	Reference Test Gain, ANSI	dB	39	—	32	—	29	—
Current	Quiescent Current	mA	1.0	1.0	1.0	1.0	1.0	1.0
	Operating Current, IEC	mA	1.0	1.0	1.0	1.0	1.0	1.0
	Operating Current, ANSI	mA	1.1	—	1.0	—	1.0	—
	Battery Type		13		312		312	
Distortion	500/800/1600 Hz typ., IEC	%	<1 / <1 / <1	<1 / <1 / <1	<1 / <1 / <1	<1 / <1 / <1	<1 / <1 / <1	<1 / <1 / <1
	500/800/1600 Hz typ., ANSI	%	<1 / <1 / <1	—	<1 / <1 / <1	—	<1 / <1 / <1	—
General Information	Frequency Range, ANSI	Hz	110-5800	—	100-7200	—	100-6000	—
	Equiv. Input Noise, IEC/ANSI <sup>1)</sup>	dB	19/18		21/20		21/20	
	Telecoil 1 mA/m 1600 Hz, IEC	dB SPL	73	82	62	71	61	70
	Telecoil HFA SPLITS, ANSI	dB SPL	97	—	90	—	83	—
Additional Information	Colors		beige / pink / brown / dark brown		beige / pink / brown / dark brown		beige / pink / brown / dark brown	
	Telecoil Boost		yes		yes		yes	
	Automatic Telecoil		optional		optional		optional	
	Push Button Modes		M, T or MT / Mute		M, T or MT / Mute		M, T or MT / Mute	
	Volume Control		yes		yes		yes	

<sup>1)</sup> Technical data measured with expansion, corresponding to Soft Noise Management™ level 3.

All measurements are made according to IEC 60118 if not otherwise mentioned. ANSI refers to ANSI S3.22-1996. The Full-On Gain setting can be programmed into the instrument from OASIS plus for verification purposes.

## Programming

Win is programmed with OASIS plus, Bernafon's outstanding fitting software. OASIS plus is a NOAH compatible, MS-Windows based PC-Fitting software. The connection is made through a HiPro or NOAHlink interface, using the appropriate programming cables and adapters.

System requirements:

NOAH 3 with HiPro or NOAHlink, programming cables and programming adapters.

# Win information

Win ITE								
			315 ITC DM		301 ITC		411 MC 401 CIC	
			IEC 60118-7 (2cc)	IEC 60118-0 (Ears.)	IEC 60118-7 (2cc)	IEC 60118-0 (Ears.)	IEC 60118-7 (2cc)	IEC 60118-0 (Ears.)
Output	OSPL 90, Peak	dB SPL	113	122	110	120	104	114
	OSPL 90, 1600 Hz	dB SPL	104	113	103	111	97	105
	HFA OSPL 90, ANSI	dB SPL	106	—	104	—	98	—
Gain	Full-On Gain, Peak	dB	43	54	40	50	34	45
	Full-On Gain, 1600 Hz	dB	34	43	31	39	25	33
	HFA Full-On Gain, ANSI	dB	36	—	32	—	27	—
	Reference Test Gain, IEC	dB	27	36	23	32	18	26
	Reference Test Gain, ANSI	dB	31	—	29	—	23	—
Current	Quiescent Current	mA	1.0	1.0	1.0	1.0	0.8	0.8
	Operating Current, IEC	mA	1.0	1.0	1.0	1.0	0.8	0.8
	Operating Current, ANSI	mA	1.0	—	1.0	—	0.8	—
	Battery Type		312		312		10	
Distortion	500/800/1600 Hz typ., IEC	%	<2 / <1 / <1	<2 / <2 / <1	<1 / <1 / <1	<1 / <1 / <1	<1 / <1 / <1	<1 / <1 / <1
	500/800/1600 Hz typ., ANSI	%	<2 / <1 / <1	—	<1 / <1 / <1	—	<1 / <1 / <1	—
General Information	Frequency Range, ANSI	Hz	100–6200	—	100–6000	—	100–6000	—
	Equiv. Input Noise, IEC/ANSI <sup>1)</sup>	dB	23/22		21/20		22/21	
	Telecoil 1 mA/m 1600 Hz, IEC	dB SPL	68	77	61	70	—	—
	Telecoil HFA SPLITS, ANSI	dB SPL	89	—	83	—	—	—
Additional Information	Colors		beige / pink / brown / dark brown		beige / pink / brown / dark brown		beige / pink / brown / dark brown	
	Telecoil Boost		yes		yes		—	
	Automatic Telecoil		optional		optional		—	
	Push Button Modes		M, DM, T or MT/ Mute		M, T or MT / Mute		M / Mute (optional)	
	Volume Control		yes		—		—	
	AI-DI	0° / max. dB	4.8/6.2		—		—	

<sup>1)</sup> Technical data measured with expansion, corresponding to Soft Noise Management™ level 3.


All measurements are made according to IEC 60118 if not otherwise mentioned. ANSI refers to ANSI S3.22-1996.

The Full-On Gain setting can be programmed into the instrument from OASIS plus for verification purposes.

Directionality measured on Kemar®. AI-DI relative to 0° and relative to maximal polar sensitivity, calculated for the horizontal plane.

# Frequency responses

**IEC 60118-7  
2cc Coupler  
(IEC 60126)**

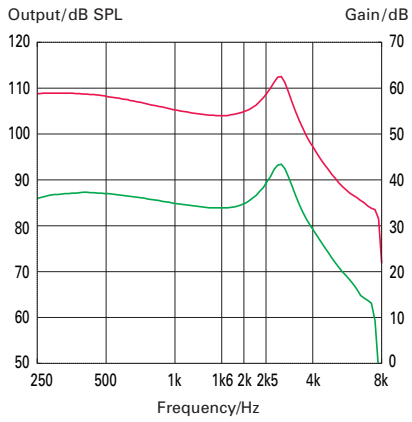


**Win 315 ITC DM**

Output OSPL 90  
Full-On Gain





**IEC 60118-0  
Earsimulator  
(IEC 60711)**

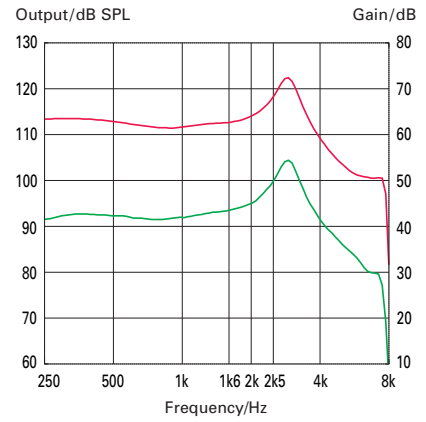



**Win 301 ITC**

Output OSPL 90  
Full-On Gain

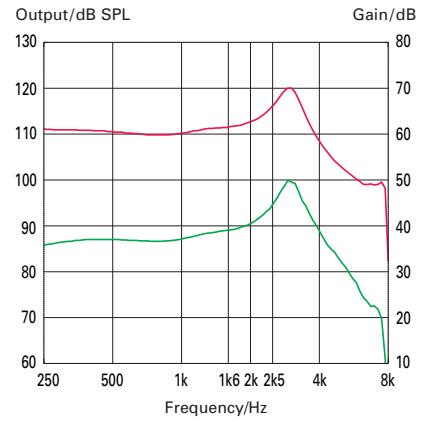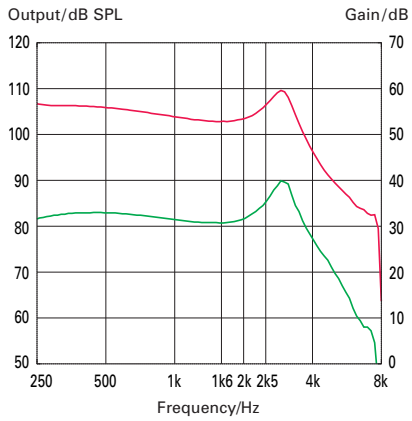

**Win 411 MC  
Win 401 CIC**


Output OSPL 90  
Full-On Gain

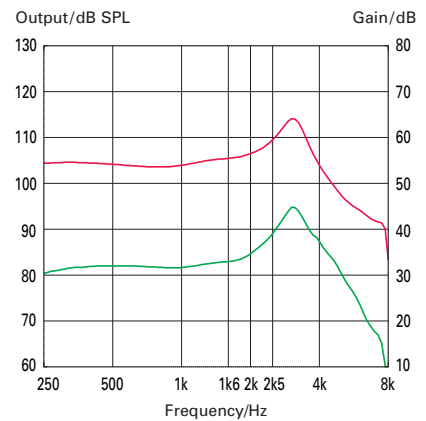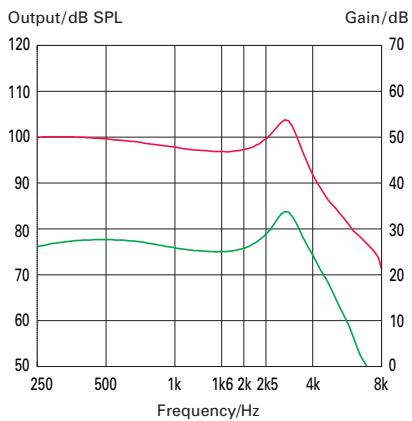



# Frequency responses


**IEC 60118-7  
2cc Coupler  
(IEC 60126)**

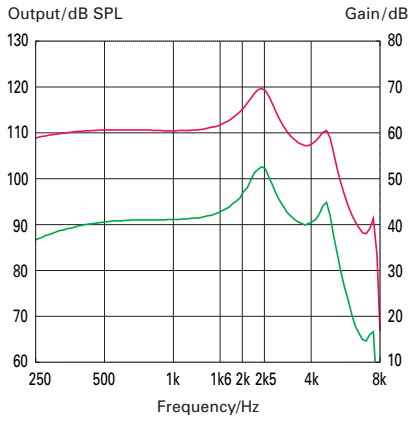

**Win 202 ITE**



Output OSPL 90  
Full-On Gain

**IEC 60118-0  
Earsimulator  
(IEC 60711)**

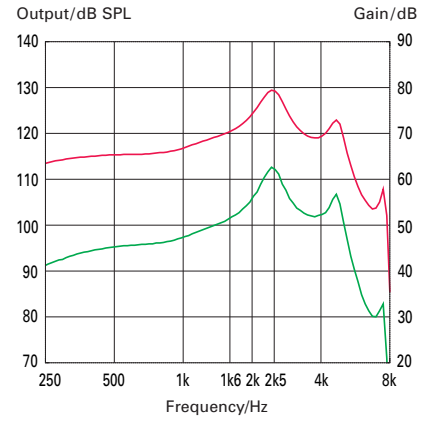



**Win 322 ITC**

Output OSPL 90  
Full-On Gain

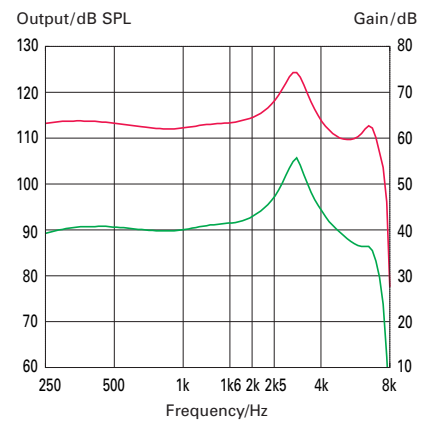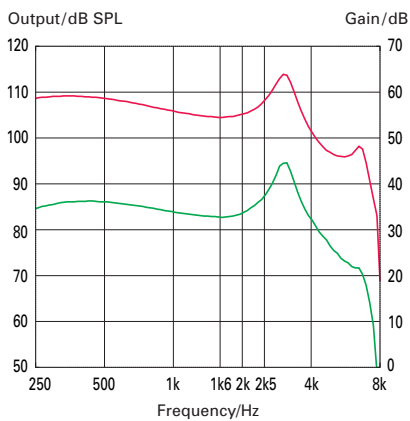

**Win 302 ITC**

Output OSPL 90  
Full-On Gain

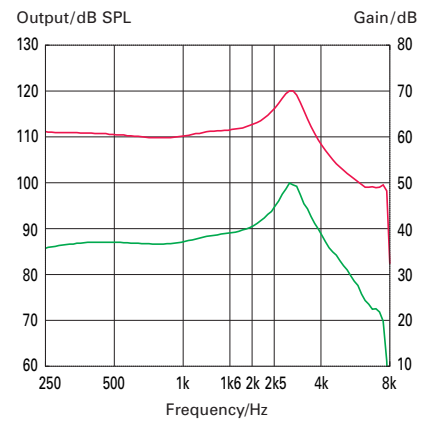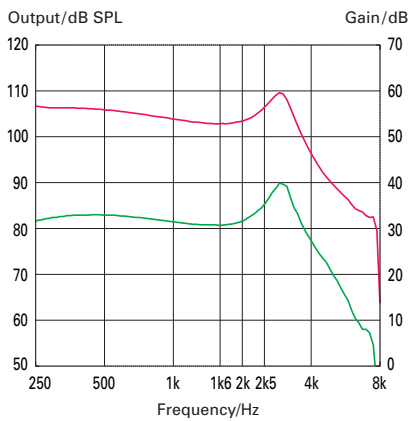

# Fitting in OASIS plus

**A high performance fitting software platform to ensure a perfect fit.**

All Win hearing instruments are fit with OASIS plus, the reliable hearing instrument fitting software from Bernafon.

Intelligent software makes your job easier.

- First fit in as few as 3 mouse clicks
- Fine control for each channel for greater adjustment precision
- 3 MPO controls
- Interactive fine-tuning to support the dialog with your client
- Efficient Feedback Management
- COSI for counseling and ensuring that your client's goals are met
- NOAHlink compatibility
- Standard cables and adapters
- Unified fitting flow with all Bernafon digital hearing instruments



**Win and OASIS plus, the perfect combination of affordable hearing instruments and high performance fitting software.**

## Programming cable attachment

### Programming cables:

Nr. 2, NEW STANDARD (HiPro)  
Blue, left Ref. 384-20-033-00  
Red, right Ref. 384-20-032-00



### Programming adapter for all ITE's

FlexConnect Ref. 390-01-180-05

### Programming cables:

Nr. 2, NEW STANDARD (NOAHlink)  
Blue, left Ref. 384-20-035-00  
Red, right Ref. 384-20-034-00


### Manufacturer:

**Bernafon AG**  
Morgenstrasse 131  
3018 Bern  
Switzerland  
Phone +41 31 998 15 15  
Fax +41 31 998 15 90  
www.bernafon.com

**Bernafon UK**  
Cadzow Industrial Estate  
Low Waters Road  
Hamilton  
ML3 7QE Scotland  
Phone +44 1698 28 59 68  
Fax +44 1698 42 14 56  
www.bernafon.co.uk

**Bernafon Australia Pty. Ltd.**  
12/97 Castlemaine Street  
Milton QLD 4064  
Australia  
Freecall 1800 809 111  
Phone +61 7 3250 0300  
Fax +61 7 3250 0372  
www.bernafon.com.au

**Bernafon New Zealand Ltd.**  
Level 1, Building F  
27-29 William Pickering Drive  
Albany, Auckland  
New Zealand  
Phone +64 9 920 42 40  
Fax +64 9 920 42 43  
www.bernafon.co.nz

**SWISS**   
Engineering

**bernafon**   
Your hearing • Our passion